

— The —  
**Fastest Laser  
Hair Removal  
Treatment on  
the Market!\***



Smart Clinic Alma  
**Soprano Titanium™**  
*Special Edition*  
English

Soprano Titanium  
*Special Edition*

\*Validated by 'Sapio Research', an independent global research agency, when compared with leading global laser hair removal brands.

# Hair Removal, Remastered. **Now, Even Faster.**



Engineered for efficiency like never before, the special edition of Soprano Titanium reduces treatment times to a fraction of what you might expect.

# 20%

speed boost\* –  
it's not just an  
improvement;  
**it's a game-  
changer.**

With Alma's revolutionary technology, the special edition of Soprano Titanium leverages the benefits of three combined wavelengths, to take the concept of speed to a whole new level. This game-changing platform enables you to treat your customers in shorter sessions, leaving them delighted with the results, while giving your practice a competitive edge.

\*Valid in SHR In-Motion mode and high fluence levels.



## Ready to **MAXimize** Results?

**Trio MAX** is a revolutionary applicator that's all about speed, efficiency, and an unparalleled customer experience. This unique applicator features a 4cm<sup>2</sup> spot size, offering an impressive increase of up to 20% in treatment speed\*.

What's more, its optimal optics delivery and enhanced efficiency in advanced sessions further elevate its performance to unprecedented levels. The Trio MAX also features a distinctive cooling system, ensuring that treatments are virtually painless.




# Synergetic Technologies. Superior Results.

## SHR

The special edition of Soprano Titanium features Alma's proprietary SHR technology, the only clinically-proven, virtually painless method of laser hair removal. Its unique combination of gradual heating, high frequency pulse delivery and strong cooling enables effective and comfortable treatments, for all skin types.

## 3D - Exclusive triple wavelength technology

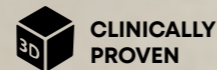

In the first solution to combine the three optimal laser wavelengths of 755nm, 810nm and 1064nm in one simultaneous emission, 3D targets different tissue depths and structures within the follicle, for maximum results.

## ICE Plus™

This advanced technology continuously cools the skin and facilitates range-controlled temperature throughout the entire treatment. Patients enjoy a virtually painless, cool and comfortable experience during the session.



## THE WAVELENGTHS

### 755 nm

Based on the Alex laser, the 755nm offers strong absorption by the chromophore melanin. This makes this the ideal wavelength in the removal of fine and light-colored hair.

### 810 nm

Good absorption in melanin and relatively deep penetration provide a balance of safety and efficacy.

### 1064 nm

Characterized by lower melanin absorption, the 1064nm wavelength is an ideal solution for darker skin types. It also offers deep penetration, targeting the bulb and papilla of the hair follicle. With water and hemoglobin absorption, the 1064nm wavelength increases the thermal profile of the overall treatment, leading to highly effective hair removal, especially for deep hair follicles.



# Manage Your Clinic. **Maximize Your Potential.**

A digital, cloud-based business tool, Smart Clinic gives you continuous access to live data and analytics, from any device. Gain valuable insights into platform performance, treatment efficiency, ROI and more, to optimize your practice.



## Additional Applicators



### TRIO 4cm<sup>2</sup>

With its large 4cm<sup>2</sup> spot size, and integration of three wavelengths, this applicator facilitates increased treatment coverage, while making the process easier and safer for patients. With TRIO 4cm<sup>2</sup>, hair removal is quick, efficient and convenient.



### TRIO - Three In One

This applicator enables the delivery of three highly effective wavelengths in one, making it suitable for all skin types. Covering a spot size of 2cm<sup>2</sup>, it simultaneously targets different tissue depths and anatomical structures within the hair follicle, increasing treatment efficiency.



### SPEED™ 810

The Speed 810 incorporates a highly effective and safe wavelength in a comfortable 2cm<sup>2</sup> spot size.



### ALEX 755

Offering powerful absorption, combined with continuous cooling, the Alex 755 enables the treatment of light and fine hair for skin types I-V.



### FACIAL TIP

Easily treat areas that are normally hard to reach, including the ears, nostrils and glabella, with the precise facial tip.



# Soprano Titanium *Special Edition*

As a leading company in the field of medical esthetics, we enable thousands of physicians to help millions of patients experience truly remarkable transformations. Alma enables you to provide the most effective, quick and safe treatment possible to each person that walks into your clinic. In a world that is constantly changing, one thing remains constant: patients need to know they are getting the best care possible – before, during and after their treatment. When our solution is in your hands, that promise becomes a reality.

Alma Lasers GmbH

Nordostpark 100-102  
90411 Nuremberg, Germany  
Tel. + 49 911 / 89 11 29-0  
Fax + 49 911 / 89 11 29-99

Email: [info@almalasers.com](mailto:info@almalasers.com)

[www.almalasers.com](http://www.almalasers.com)

© 2025 Alma. All rights reserved.  
Alma, its logo, Soprano Titanium, TRIO, Trio MAX, ICE Plus, Smart Clinic, SHR & Speed are trademarks and/or registered trademarks of Alma in the United States and/or other countries.

PBST19112301\_03

CONNECT WITH ALMA

